

Vortex Mixer

VT-02 I VT-03

Usage - comfortable everyday mixing operations

Maximum Speed - adjustable up to 2800 rpm (7-110-00-0) / 4200 rpm (7-120-00-0)

Comfortable use - selection of different modes (touch, stand-by and continuous), depending on the desired use / smooth and silent operation

Stable operation - through inbuilt counter balance system

Digital display (VT-03) – adjustable speed setting with interval of 10rpm and time setting from 0-999 mins / display speed and time alternatively every 5 seconds

Programmable pulse mode – On time and OFF time can be individually modified


Multiple applications - optional attachments for micro plate, micro tubes and flask sticker for permissible speed for various attachments

Robust design – brushless DC motor enables maintenance-free and long-lasting operation


Small footprint saves valuable bench space - for use in hoods in the temperature range from 5 to 40°C



Optional attachments*

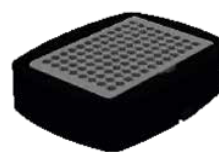

Foam for micro plate

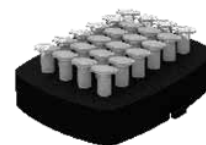

Foam for microtubes



Rubber holder for flask

Ordering information

Description	Cat. No.	Price w/o VAT
AHN myLab® VT-02 Vortex Mixer 2800 rpm	7-110-00-0	195.00 €
AHN myLab® VT-03 Digital Vortex Mixer 4200 rpm	7-120-00-0	280.00 €
Universal platform	7-110-01-0	20.00 €
Foam for microtubes*	7-110-02-0	22.00 €
Foam for micro plate*	7-110-03-0	27.00 €
Rubber holder for flask*	7-110-04-0	22.00 €

*only applicable in combination with the universal platform 7-110-01-0

Specifications	
Orbital diameter	4 mm
Speed accuracy	± 5.0%
Max. loading capacity	500 g
Protection class	IP21
Power supply	220 240 VAC, 50 Hz
Dimensions (WxDxH)	205x136x138.5 mm
Weight	3.0 kg

A&E Scientific srl – bv

Rue de Lekernay, 11 / Lekkernijstraat, 11
B-7850 Marcq / Mark

Tel : +32(0)2 396 19 10 Fax : +32(0)2 396 06 73
info@aescientific.be www.aescientific.be